

iCare IC100



Easy to use with minimal training

iCare
For better perception

iCare IC100 tonometer

The iCare's patented rebound technology is a reliable technique for obtaining accurate and rapid IOP measurements from all your patients. Measuring requires no anesthetic drops, air puffs or other preparation.

No anesthesia, no air puff

The easy-to-use iCare IC100 tonometer incorporates guiding features and an intuitive user interface to maximize efficiency in clinical routines. No specialized skills are needed. Just load, align, and measure. The accurate iCare IC100 tonometer is the professional's choice.




Key features and benefits



iCare AMS
Dual-mode operation (single or rapid 6-series measurement)



iCare EasyNav
Advanced navigation interface

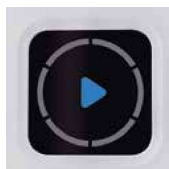

OLED display
A bright OLED display screen guides operation and presents results



iCare EasyPos
Alignment indicator: green confirms tonometer positioning



Functionality with comfort

-  Effortless probe loading using the original iCare probe.
-  The momentary contact of the probe with the cornea is harmless and barely noticed by the patient.
-  A probe cannot accidentally drop from the device.
-  With iCare IC100 IOP measuring is fast, accurate and reliably repeatable.
-  Simple user menu, available in multiple languages.
-  All six measurements required for calculation of a result can be performed individually or in an automated series, at the press of a button.

Technical data

iCare IC100	
Model name	TA011
Dimensions (W x H x D)	29 mm x 95 mm x 215 mm / 1.1" x 3.7" x 8.5"
Weight	135 g (without batteries) / 0.3 lbs (without batteries) 230 g (with 4 x AA batteries) / 0.5 lbs (with 4 x AA batteries)
Power supply	4 x AA non-rechargeable batteries, 1.5 V alkaline LR6
Measurement range	7 mmHg – 50 mmHg
Accuracy	± 1.2 mmHg (≤ 20 mmHg) and ± 2.2 mmHg (> 20 mmHg)
Repeatability	< 8 % (coefficient of variation)

iCare. For better perception.

iCare is a trusted partner in ophthalmic diagnostics, offering physicians fast, easy-to-use, and reliable tools for diagnosis of glaucoma, diabetic retinopathy, and macular degeneration (AMD). Our product assortment includes automated TrueColor imaging devices, perimeters and handheld rebound tonometers.

We believe that ophthalmic care should be accessible, effortless, and reliable, and we aim to establish the next level of eye care.



iCare Finland Oy

Äyritie 22
01510 Vantaa, Finland
Ph. +358 9 8775 1150
info@icare-world.com

Centervue S.p.A.

Via San Marco 9H
35129 Padova, Italy
Ph. +39 049 501 8399
info@icare-world.com

iCare USA, Inc.

4700 Falls of Neuse Rd. Ste 245
Raleigh, NC. 27609
Ph. +1 888.422.7313
Fax +1 877.477.5485
infoUSA@icare-world.com

www.icare-world.com

icare

For better perception

iCare is a registered trademark of Icare Finland Oy. Centervue S.p.A., Icare Finland Oy and Icare USA, Inc. are parts of Revenio Group and represent the brand iCare. Not all products, services or offers referenced in this brochure are approved or offered in every market and approved labeling and instructions may vary from one country to another. Product specifications are subject to change in design and scope of delivery and as a result of ongoing technical developments.